Gewa IR-1SP Micro

SE: IR-sändare

GB: IR transmitter

NO: IR-sender

FR: IR-Émetteur

DE: IR-Sender

Art.no 401104

ABILIA

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GB: English Manual Gewa IR-1SP Micro

Introduction

The Gewa IR-1SP Micro is a programmable IR transmitter, which can learn one channel from another remote control for TV, DVD, stereo etc.


The Gewa IR-1SP Micro is a part of the “GewaLink-system” and therefore can be used together with all of these products.

General information about Gewa IR-1SP Micro

Before using your IR-1SP Micro there are certain things that you should know:

- On delivery the GewaLink channel 16 is factory programmed, meaning that it is possible to use the transmitter immediately after the batteries has been inserted. See section *Inserting Batteries*.
- For programming press the programming button **P** (using a pen tip or other narrow object) at the **same time** as you press the function button **O**.
- If you press the wrong button or make a mistake when programming you can always press the programming button **P**. It works as a *Reset*-button.
- It is important not to wait too long between the different steps when programming. After a certain time your IR-1SP Micro automatically leaves the programming mode to save power.
- ^{3.0} as a superscript indicates the IR-1SP Micro version.

This picture shows the various terms and functions which are described in the text.

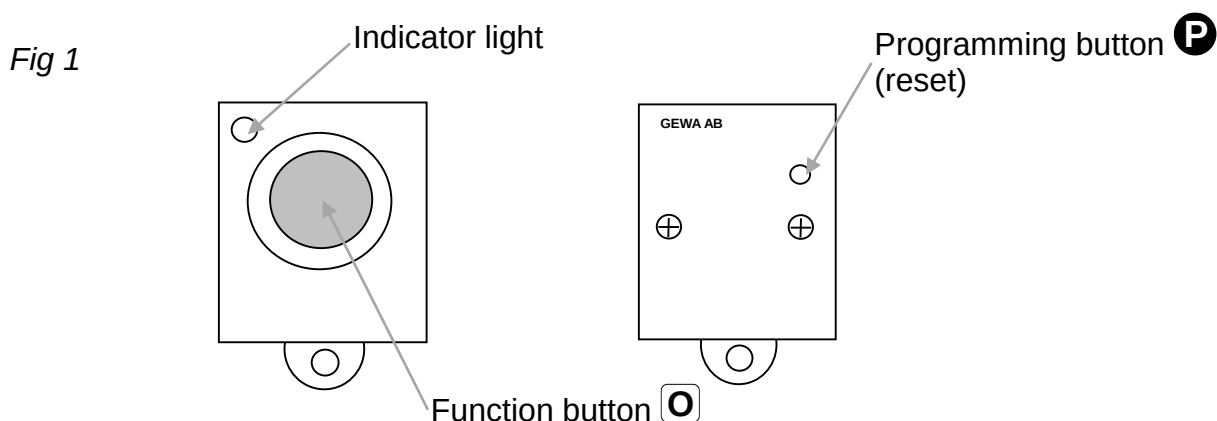The following symbols appear in the text throughout this manual:

= programming button

= function button on the IR-1SP Micro


= Any button on another transmitter to be programmed

Inserting Batteries

Unscrew the screws holding the bottom of the IR-1SP Micro using a Phillips screwdriver. Remove the bottom and insert three 1.5V silver oxide batteries type SR44, or 401930 Battery pack for Micro and Wrist transmitters, available from Abilia or Abilia dealer. Insert the batteries according to the guidelines at the bottom of the battery compartment. Pay attention to the polarity! It is important to use specified batteries for maximum transmission range.

When you have inserted the new batteries, **press the programming button **P**** on the underside.

Fig 2

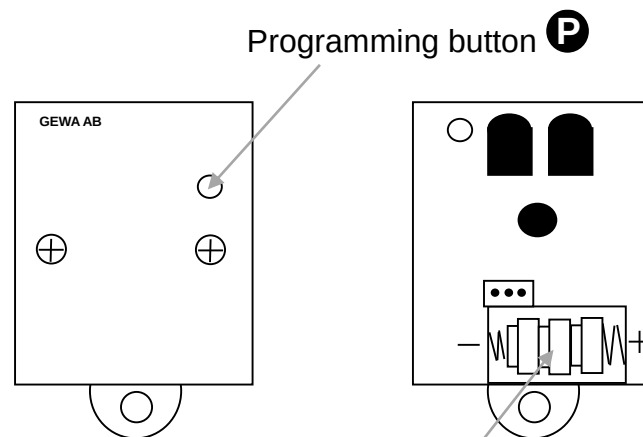Batteries: 3 x SR44
or 401930 Battery pack for Micro and Wrist transmitters

Recording IR Channel from another Transmitter


Make sure the batteries in both transmitters are fully charged. Do not use the transmitters in intense light (e.g., sunlight) when programming.

Place the two transmitters on a table with the original transmitter directed towards the front of your IR-1SP Micro, see picture.

1. Press the programming button **P** and the function button **O** at the same time.

To check that the transmitters are in the right position, keep any button pressed on the original transmitter and move the original transmitter to the left and right until the indicator light goes out. Then move the original transmitter backwards until the indicator light goes out (one metre for some transmitters). Place the original transmitter at the centre of these maximum range positions. Certain types of transmitters only emit a short flash each time the button is pressed. In these cases, press several times.

2. Press and release the function button **O** on the IR-1SP Micro. The indicator light now emits a steady red light.
3. Press and hold the button **Y** on the original transmitter that shall be recorded to IR-1SP Micro. Release the button when the indicator light goes out. It might take up to 10 seconds. If the indicator light turns red again then press the button again until it goes out.
4. Now test your recorded function on the device that is to be controlled.

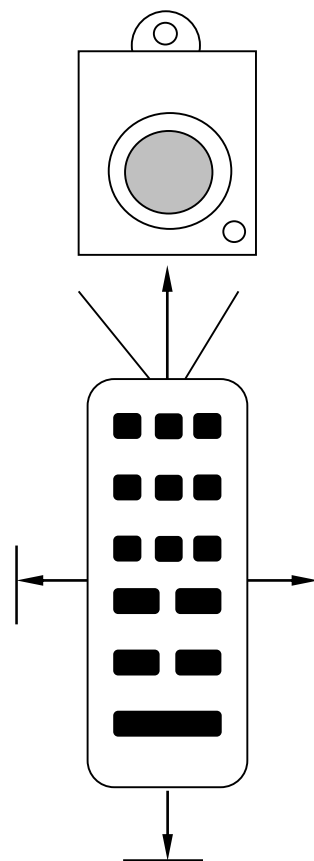Recording IR Channel

Press **P** + **O**, **O**, select button **Y**

^{3.0} Tip! IR-1SP Micro can record codes from most IR transmitters on the market with a few exceptions. You may have problems with a few or all of the functions for a certain IR transmitter. In these cases it is possible to use another recording method. Instead, do as follows in step 1; Press the programming button and the function button at the same time but keep them pressed for approximately 2 seconds. IR-1SP Micro is now in another recording mode. Otherwise use the same programming procedure.

Mounting Gewa IR-1SP Micro

The IR-1SP Micro can be used hanging on a strap around the neck. Thread the thin loop of the neck strap through the attachment loop on the transmitter. Adjust the length of the neck strap.

Warning!

The string lock is a safety feature to prevent strangling. Make sure that the ends of the neck strap can be pulled out eg. if the transmitter gets caught on something. Never make knots in the neck strap!

Cleaning

Your IR-1SP Micro can be cleaned with a damp cotton wool pad with surgical spirit, Chlorhexidine or ethyl alcohol.

It can also be cleaned with a damp cloth (weak solution of washing up liquid). Wipe dry after cleaning.

Warning!

Do not use any disinfectant except surgical spirit, Chlorhexidine or ethyl alcohol.

These liquids are inflammable, handle with care!

Operating Environment

The IR-1SP Micro is designed for use indoors and outdoors but should be protected against splashing liquid.

Technical data

Technical data

Power consumption at rest:	Approx. 10 μ A
Power consumption at transmitting:	Approx. 10 mA (GewaLink channel)
Battery life:	Approx. 1 year (30 seconds transmitting per day with GewaLink channel)
Operating force button:	Approx. 170g
Range:	10-15m (GewaLink channels)
Measurements (L x W x H):	51 x 36 x 16 mm
Weight (incl. batteries):	30g
Operating temperature (excl. batteries):	-20°C to +45°C
Product life:	6 years

Accessories included with your Gewa IR-1SP Micro

Art.no	Name	Description	Qty
890410	Neck strap	Allows IR-1SP Micro to be carried around the neck	1 pcs
419300	Battery pack for Micro and Wrist transmitters	Battery pack 1,5 V Silver oxide	1 pcs

Safety and responsibility

Abilia has conducted a risk analysis for the product with reference to SS-EN 14971 for Medical equipment. The instructions for use and the technical specifications for the product have been drawn up in respect of the safety assessments carried out in this document.

It is always important to be aware of safety when first using this product and with its subsequent use. It is not possible to describe in the instructions all the situations in which the product must not be used. If you are not sure whether you have started using the product properly or set it up correctly, please contact your supplier.

Abilia is not liable for any loss, damage or consequential loss resulting from any use, correct or otherwise, lack of follow-up, incorrect set-up, incorrect connection or similar misuse. Abilia financial liability is at all times limited to the value of the product itself.

Recycling Waste Electrical and Electronic Equipment

The IR-1SP Micro contains electronics components that must be recycled. When the product is no longer of use, hand it in to a recycling centre for waste electrical and electronic equipment.

Conformity with the Medical Devices Directive

The IR-1SP Micro is marketed as a technical aid for people with disabilities and meets the requirements in the Medical Devices Directive 93/42/EEC including relevant EMC and electrical safety requirements.

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Abilia AB, Kung Hans väg 3, Sollentuna, Box 92, 191 22 Sollentuna, Sweden